

Lot Conversion

Lot Nxxxxxxx PT

ACL TOP 770s - 23080002

03/28/2024

Prepared for

Inspection

Prepared by
me

Authorized Signer Information:

Signed by: ppatel2

Signed date: 03/28/24

Authorized Signer Information

Signed by: ppatel2

Signed date: 03/28/24


Summary

Active Lot	Test	N	Slope	Intercept	Status
N0239577	PT RecombiPlasTin	21/20	1.326	-5.308	Pass

Authorized Signer Information

Signed by: ppatel2

Signed date: 03/28/24

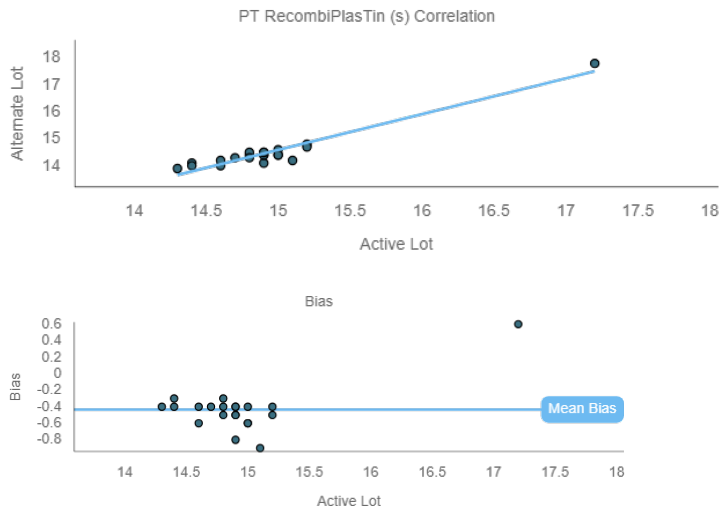

Parameter	Ac.Lot: N0239577	Al.Lot: N0824895
N	21	21
Mean	14.93	14.50
SD	0.58	0.79
2SD	1.15	1.58
Median	14.90	14.40
Geometric Mean	14.92	14.48
Low Range	13.78	12.92
High Range	16.08	16.08

Parameter	Calculated	Acceptance Criteria
N	21	20
R	0.965	
Slope	1.326	
Intercept	-5.31	
Mean N0239577	14.93	
Mean N0824895	14.50	
Mean of Means	14.72	
Mean of Bias	-0.43	
Mean of Bias %	-2.94	
SD of Bias	0.28	
Mean % Difference	-2.95	
Allowable Limit	-2.90	-

Authorized Signer Information

Signed by: ppatel2

Signed date: 03/28/24

Status	Pass	
--------	------	--

Authorized Signer Information

Signed by: ppatel2

Signed date: 03/28/24

Results

#	Container	N0239577	N0824895	Mean	Raw Diff.	% Diff.
1	1	14.90	14.10	14.50	-0.80	-5.37
2	2	14.60	14.00	14.30	-0.60	-4.11
3	3	14.40	14.10	14.25	-0.30	-2.08
4	4	14.70	14.30	14.50	-0.40	-2.72
5	5	15.00	14.40	14.70	-0.60	-4.00
6	6	15.10	14.20	14.65	-0.90	-5.96
7	7	14.80	14.40	14.60	-0.40	-2.70
8	8	15.20	14.80	15.00	-0.40	-2.63
9	9	14.90	14.40	14.65	-0.50	-3.36
10	10	14.90	14.40	14.65	-0.50	-3.36
11	11	14.40	14.00	14.20	-0.40	-2.78
12	12	14.60	14.20	14.40	-0.40	-2.74
13	13	14.30	13.90	14.10	-0.40	-2.80
14	14	15.20	14.70	14.95	-0.50	-3.29
15	15	15.00	14.60	14.80	-0.40	-2.67
16	16	14.90	14.50	14.70	-0.40	-2.68
17	17	15.00	14.40	14.70	-0.60	-4.00
18	18	14.90	14.50	14.70	-0.40	-2.68
19	19	14.80	14.50	14.65	-0.30	-2.03
20	20	14.80	14.30	14.55	-0.50	-3.38
21	21	17.20	17.80	17.50	0.60	3.49

Authorized Signer Information

Signed by: ppatel2

Signed date: 03/28/24