

Reference Interval/NPP

ACL TOP 570 - 2345678

ACL TOP 770 LAS - 1234567

07/30/2024

Prepared for

Prepared by

Authorized Signer Information:

Signed by: admin

Signed date: 07/30/24

Authorized Signer Information

Signed by: admin

Signed date: 07/30/24



Summary

Lot	Test	N Results	Reference Range	Status
N0239577	PT RecombiPlasTin	40/10	14.33 to 15.36	Pass
N0824895	PT RecombiPlasTin	40/10	13.13 to 14.83	Pass

Authorized Signer Information

Signed by: admin

Signed date: 07/30/24

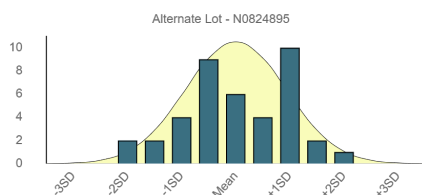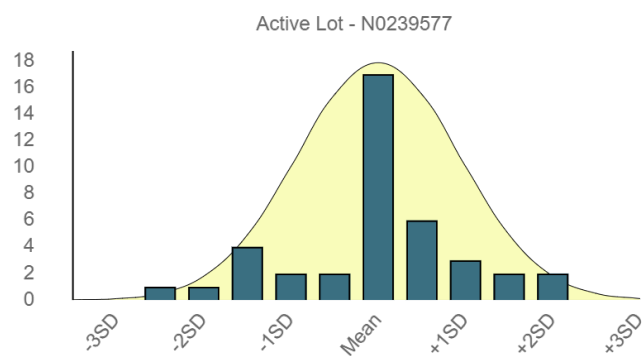

Parameter	Ac.Lot: N0239577	Al.Lot: N0824895	Acceptance Criteria
N	40	40	10
Mean	14.84	13.98	
SD	0.26	0.42	
2SD	0.51	0.85	
Median	14.90	13.95	
Geometric Mean (NPP)	14.84	13.97	
Geometric SD	1.02	1.03	
Low Range	14.33	13.13	
High Range	15.36	14.83	
Values in Reference Interval	0.00	0.00	-
Status	Pass		

Authorized Signer Information

Signed by: admin

Signed date: 07/30/24

Results

#	Container	N0239577	Instrument	N0824895	Instrument
1	1	14.90	ACL TOP 570	14.10	ACL TOP 570
2	2	14.60	ACL TOP 570	14.00	ACL TOP 570
3	3	14.40	ACL TOP 570	14.10	ACL TOP 570
4	4	14.70	ACL TOP 570	14.30	ACL TOP 570
5	5	15.00	ACL TOP 570	14.40	ACL TOP 570
6	6	15.10	ACL TOP 570	14.20	ACL TOP 570
7	7	14.80	ACL TOP 570	14.40	ACL TOP 570
8	8	15.20	ACL TOP 570	14.80	ACL TOP 570
9	9	14.90	ACL TOP 570	14.40	ACL TOP 570
10	10	14.90	ACL TOP 570	14.40	ACL TOP 570
11	11	14.40	ACL TOP 770 LAS	14.00	ACL TOP 770 LAS
12	12	14.60	ACL TOP 770 LAS	14.20	ACL TOP 770 LAS
13	13	14.30	ACL TOP 770 LAS	13.90	ACL TOP 770 LAS
14	14	15.20	ACL TOP 770 LAS	14.70	ACL TOP 770 LAS
15	15	15.00	ACL TOP 770 LAS	14.60	ACL TOP 770 LAS
16	16	14.90	ACL TOP 770 LAS	14.50	ACL TOP 770 LAS
17	17	15.00	ACL TOP 770 LAS	14.40	ACL TOP 770 LAS
18	18	14.90	ACL TOP 770 LAS	14.50	ACL TOP 770 LAS
19	19	14.80	ACL TOP 770 LAS	14.50	ACL TOP 770 LAS
20	20	14.80	ACL TOP 770 LAS	14.30	ACL TOP 770 LAS
21	100	14.80	ACL TOP 570	13.20	ACL TOP 570
22	101	14.90	ACL TOP 570	13.70	ACL TOP 570
23	102	14.50	ACL TOP 570	13.50	ACL TOP 570
24	103	14.80	ACL TOP 570	13.70	ACL TOP 570
25	104	15.00	ACL TOP 570	13.70	ACL TOP 570
26	105	14.80	ACL TOP 570	13.70	ACL TOP 570
27	106	14.80	ACL TOP 570	13.50	ACL TOP 570
28	107	14.90	ACL TOP 570	13.60	ACL TOP 570
29	108	14.90	ACL TOP 570	13.70	ACL TOP 570
30	109	15.00	ACL TOP 570	13.40	ACL TOP 570
31	200	14.20	ACL TOP 570	13.20	ACL TOP 570
32	201	14.50	ACL TOP 570	13.50	ACL TOP 570
33	202	14.70	ACL TOP 570	13.70	ACL TOP 570
34	203	15.10	ACL TOP 570	14.00	ACL TOP 570
35	204	15.00	ACL TOP 570	13.80	ACL TOP 570
36	205	14.90	ACL TOP 570	13.70	ACL TOP 570
37	206	15.10	ACL TOP 570	13.70	ACL TOP 570

Authorized Signer Information

Signed by: admin

Signed date: 07/30/24

38	207	15.30	ACL TOP 570	13.90	ACL TOP 570
39	208	15.30	ACL TOP 570	13.90	ACL TOP 570
40	209	14.90	ACL TOP 570	13.40	ACL TOP 570

Authorized Signer Information

Signed by: admin

Signed date: 07/30/24